

a_w -measurement simple as never before!

For more than 50 years Novasina is the leader in the production of precision instruments for water activity measurement in foodstuff, processed food as well as pharmaceutical, chemical products and

cosmetics. Thanks to a longstanding experience and research and development we are able to offer today reliable, high quality "Swiss made" products that set new standards in the water activity measurement.

LabTouch-aw

Smart water activity measurement device



Let's welcome the next generation of high quality **water activity meters** made by Novasina. Touch a new world, driven by the ease of use and fast menu access by a **touch-screen** combined with a **semi-temperature-stabilised** measurement chamber for consistent and reproducible measurement results (a_w -value). The latest sensor and measurement technology coupled with re-useable humidity calibration standards provide a very economical solution to precision water activity measurements. Furthermore, the standards are an easy method for double-checking the instrument and sensor's high performance.

The **LabTouch-aw**, with its ergonomic design and the touch-screen, makes water activity measurements easier than ever before. Semi-temperature-control and active sample temperature measurements, by an integrated infrared sensor, provide the highest possible reliability of measurements.

All measurement data, including the important protocols are stored on a SD card and can be transferred to a PC or printer. Special software is available for viewing results and analysis of the data. This assures confidence in readings for quality assurance and traceability of all measurement values.

In summary, the **LabTouch-aw** assures high measurement accuracy and reproducibility, long-term stability, ease of use and an economical measurement system. Get ready for the new standard for **aw-value** measurement

Advantages :

- Easy accessible menus by large touch-screen
- Semi-temperature-stabilized measurement chamber
- Re-usable SAL-T humidity standards and long-life sensor for economy
- Fast and precise aw-measurement thanks to the worldwide unique "Novalyte" measurement technology
- Specific chemical sensor protection filters available
- Calibration data is stored on the sensor
- Data logging function with SD card
- Factory calibration at 7 a_w -value points
- Checking, testing and adjusting possibilities (SAL-T humidity standards)
- Visualization and analysis software "Novalog MC"



Accurate a_w -measurement instrument
LabTouch-aw



novasina
The Art of Precision Measurement

Neuheimstrasse 12
CH-8853 Lachen SZ
Tel: +41 (0)55 642 67 67
Fax: +41 (0)55 642 67 70
info@novasina.ch / www.novasina.com

Specifications LabTouch-aw :

	A_w-value	Sample temperature
Measuring principle :	Resistive - electrolytic	Surface Infrared-measurement + NTC
Measuring range :	0.03 1.00 a _w	5...45 °C (41.....113°F)
Accuracy:	+/- 0.005 a _w	+/- 0.15 °K
Range :	0.10.... 0.97 a _w (10 ...97% rh)	0... +50 °C
Resolution :	0.001 a _w	0.1 °C
Accuracy of temp. control :	+/- 0.05°C	Semi-temperature control, programmed temperature must be always +2°C above environmental temperature!
Mains supply :	5 VDC +/- 6%	max. power requirement : 10 W
Power supply :	90....264 VAC	50 / 60 Hz, output -> +5 VDC
Display :	4.3"	color LC-Display
Operating :	System start	-> open chamber Touch screen and easy accessible menu structure
Communication :	SD card	type : SD / SD _{HC} data system: FAT-16 / FAT-32
Housing :		Two-part PC/ABS design housing measurement chamber aluminium anodized
Weight:	2.8 kg	
Protection class :	IP 22	Dimension approx. : 300 x 200 x 105 mm
Volume	21.1 ml	standardised sample dishes
Measuring chamber	40 x 12 mm	Spring-loaded measurement head
Special feature :	Intelligent „Novalyte“ measurement sensor CM-2 with calibration point memory	






Accessories + Options

	Article-Nr.
SD memory card	260 0193
Novalog MC visualisation for PC	260 0188
Replacement sensor CM-2	111 9983
Power supply universal	260 0505
Protection filter REDOX	111 7212
Protection filter eVC-21	111 1001
Protection filter eVC-26	111 1003
Protection filter eVALC	111 0995
Pre-filter white 5 pieces	260 0183
Tension ring	260 0210
ePW sample dishes 100 pieces	111 0601

Humidity/aw Standards

	Article-Nr.
Humidity Check SAL-T 11% RH	111 0930
Humidity Check SAL-T 33% RH	111 0932
Humidity Check SAL-T 58% RH	260 0189
Humidity Check SAL-T 75% RH	111 0936
Humidity Check SAL-T 84% RH	111 9886
Humidity Check SAL-T 97% RH	111 9987
with UKAS certificate :	
Humidity Check SAL-T/C 11% RH	111 1056
Humidity Check SAL-T/C 33% RH	111 1057
Humidity Check SAL-T/C 58% RH	260 0190
Humidity Check SAL-T/C 75% RH	111 1063
Humidity Check SAL-T/C 84% RH	260 0276
Humidity Check SAL-T/C 97% RH	111 9981

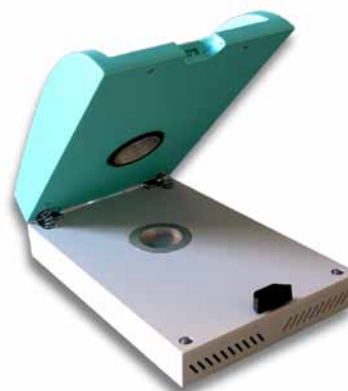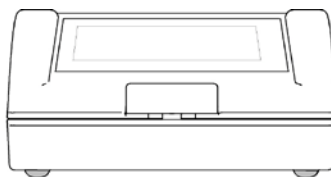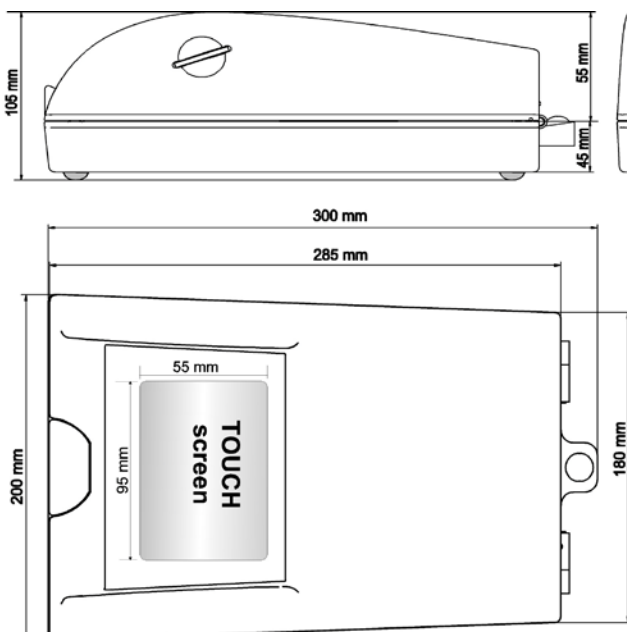

004196.00